

Simultaneous Multi-gas Analyzer Raman Spectroscopy for Fluorine Chemical Industry

Our Product Introduction

for more products please visit us on spectralanalyser.com

Basic Information

- Place of Origin: CHINA
- Brand Name: JINSP
- Certification: CE ISO9001
- Model Number: RS2600
- Minimum Order Quantity: 1
- Price: Negotiable
- Packaging Details: International Shipping Package
- Delivery Time: 90-120 working days
- Payment Terms: T/T, Western Union
- Supply Ability: 50 Unit 90-120 days



Product Specification

- Model Name: RS2600 Multi-gas Analyzer
- Laser Excitation Wavelength: $532 \pm 0.5 \text{ Nm}$
- Spectral Coverage: $200 \sim 4200 \text{ Cm}^{-1}$
- Analysis Time: 2 Seconds
- Detection Limit: Ppm Level
- Measurement Time: Real-time
- Highlight: **Simultaneous Raman Gas Analyzer, Fluorine Chemical Industry Raman Gas Analyzer, Multi-gas Raman Gas Analyzer**

Product Description

RS2600 Raman Based Multi-gas Analyzer

Engineered for aggressive environments, the RS2600 analyzer breaks barriers in fluorine chemical monitoring. Capable of detecting corrosive agents like F_2 , BF_3 , and HF with ppm-level sensitivity, this system eliminates traditional pressure control requirements while delivering instantaneous multi-gas readings. Its non-invasive design preserves sample integrity during continuous analysis of reactive compounds.

From PF_5 synthesis to HCl scrubbing systems, the platform ensures worker safety through real-time alerts. The unique wide-range detection (1ppm-100%) coupled with sub-second response makes it indispensable for process control in fluoropolymer manufacturing and semiconductor etching applications.

Specifications

Components	N_2 , H_2 , O_2 , CO , CO_2 , H_2S , $CmHn$, etc.
Analysis time	2 seconds
Measurement uncertainty	$\leq 0.2\%$
Laser excitation wavelength	532 ± 0.5 nm
Spectral coverage	200 ~ 4200 cm^{-1}
Spectral resolution	≤ 8 cm^{-1} at full spectral range
Air circuit interface	6 mm standard tube fitting (3 mm, 1/8", and 1/4" are optional)
Input voltage	100~240 VAC, 50~60 Hz
Sample gas temperature	-50 ~ 40 $^{\circ}C$
Sample gas pressure	1.0 MPa
Unit dimensions	485 mm (Width) × 350 mm (Height) × 600 mm (Depth)
Weight	40 kg

Technical highlights

Non-destructive gas analysis: Analytes including homonuclear diatomic gases (F_2 , Cl_2 , etc.) and isotopic gases (H_2 , D_2 , T_2 , etc.)

Short detection time: Data acquired in seconds


Minimal maintenance: Resistance to high pressure, direct detection without consumables (chromatographic column or carrier gas)

Wide detection range: LOD at ppm level, detection range up to 100%

Application: 1. Monitoring of F_2 , N_2 , HF and other components in the fluorination process.

Monitoring of F_2 , N_2 , HF and other components in the fluorination process


component	content
F_2	2%
N_2	79%
HF	0.2%
CO_2	3%
CO	0.0%
O_2	3%
CF_4	10%

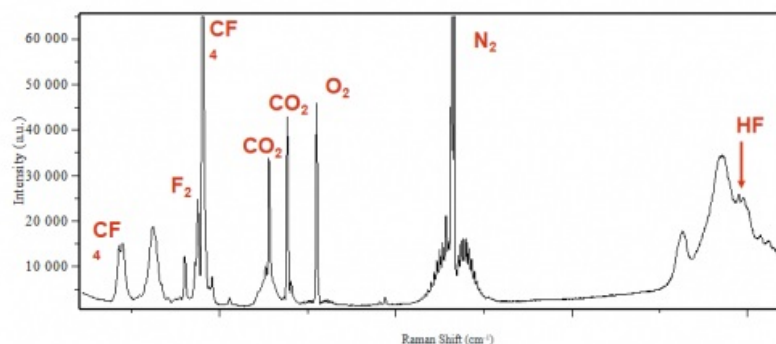


Application: 2. Quantitative analysis of different Elements in Petrochemical Natural Gas

Effective identification and quantification of various components in natural gas mixtures within 2 seconds, such as alkanes, alkenes, alkynes, etc.

Quantitative analysis of different Elements in Petrochemical Natural Gas

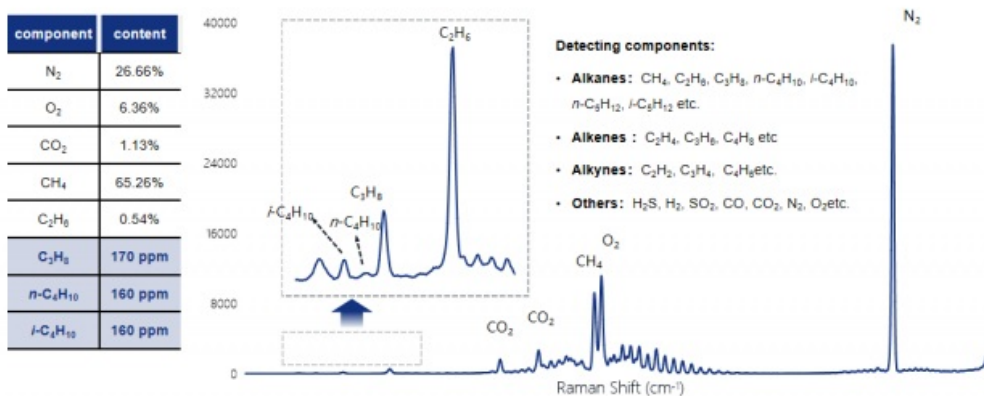


Quantitative analysis of different Elements in Petrochemical Natural Gas

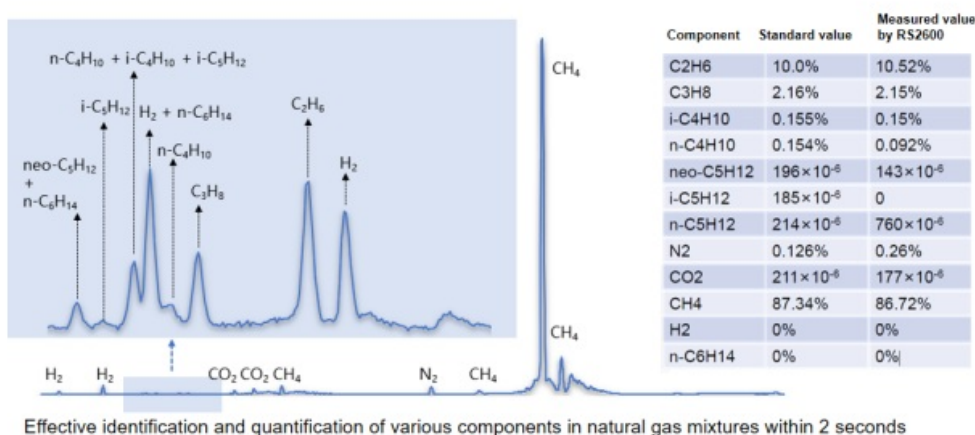


Application: 3. Detection of components in reaction and emitted gases of fluorine materials, such as reactant gases, common gases, fluoinated gases.

Detection of components in reaction and emitted gases of fluorine materials

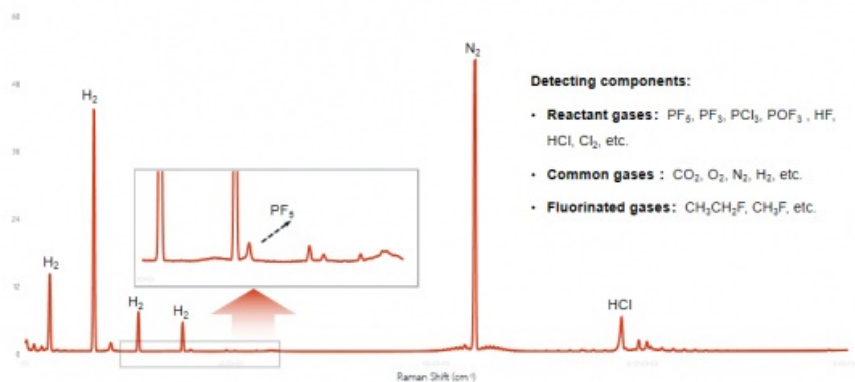


Application: 4. Detection of impurity gases and electronic special gases.

Detection of Impurity Gases in Electronic Special Gases

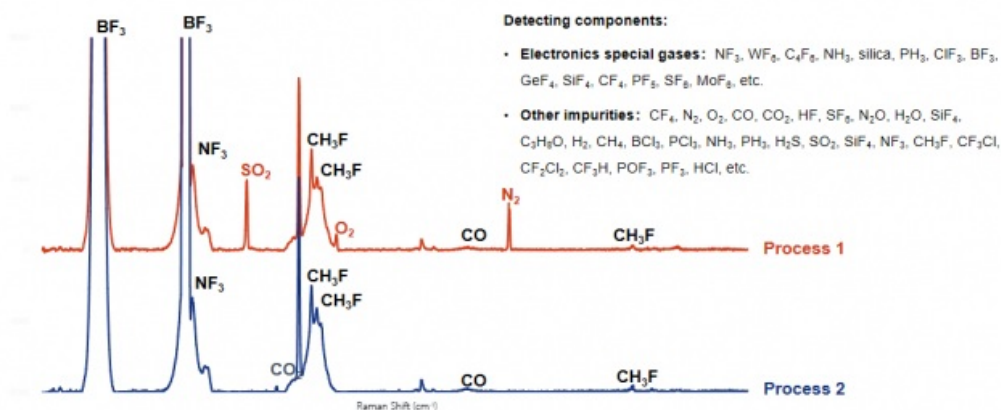

Gas Analyzer Features:

1. By constructing quantitative models of multiple standard curves and combining them with chemometric methods, the correlation between spectral signals and the content of various component substances can be effectively established. This comprehensive method not only improves the accuracy of analysis, but also enables the simultaneous processing and analysis of multiple different chemical components, providing a powerful tool for quantitative analysis of complex samples.
2. Even in the case of changes in sample gas pressure or fluctuations in testing conditions, this quantitative analysis method can still maintain the accuracy of the results. This is due to its high stability and robustness, ensuring consistent analysis results in any environment, providing reliable data support for scientific research and industrial applications.

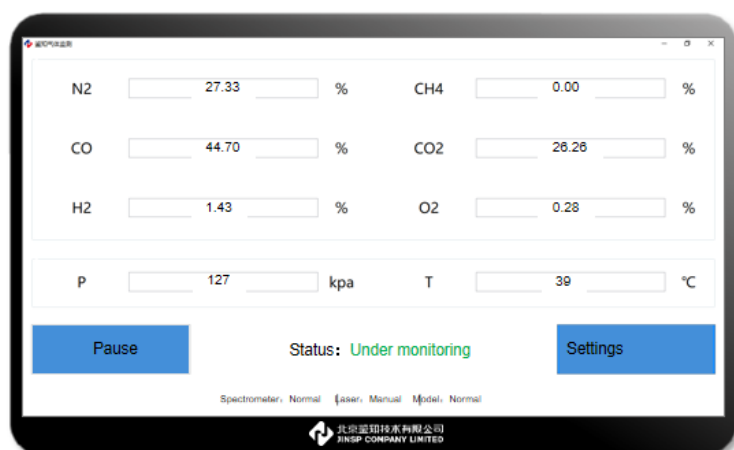

Company Introduction:

JINSP Company Limited, abbreviated as "JINSP", is a professional supplier with over 17 years of experience in spectral detection technology products, including Raman, FT-IR, LIBS technologies, etc. After 17 years of technology accumulation, the company's core key technologies have reached the international leading position at the level, and the cumulative number of patent applications exceeded 200.

JINSP offers over twenty spectroscopic products across various fields, including pharmaceutical and chemical industries, public security, and customs. Products are available nationwide and are exported to over 30 countries, with cumulative sales exceeding 3,000 units.

Benefit from 30+ R&D engineers, including 4 Ph.D., JINSP is deeply rooted in the field of personalized product customization, and is committed to meeting the diverse and unique needs of customers with excellent professional technology and innovative design capabilities.

Company Profile



Exhibition



Certifications



FAQ:

Q: This is the first time I use, is it easy to operate?

A: We will send you manual and guide video in English, it can teach you how to operate the spectrometer. Also our technicians will offer professional technical operation meetings.

Q2: Can you offer a operation training?

A: Your technicians can come to our factory for a training. Jinsp technical engineers can go to your place for local support. (installation, training, debugging, maintenance)

Q3: How to receive a best price in the shortest time?

A: When you send us an inquiry, please kindly offer details with wavelength, detector, effective pixels, focal length and so on. We will send you quotation with details soon to your email.

Q: Do you have MOQ requirement of each order?

A: No, we support order for even 1pc order.

Q: What is the product warranty?

A: All product has 1 year warranty, support extending warranty period with extra cost. Please contact with our sales person.



JINSP Company Ltd.



8618620854039



phoebeyu@jinsptech.com



spectralanalyser.com

21st Floor, Building D, Tsinghua Tongfang Science and Technology Plaza, Haidian District, Beijing China