

200nm-1100nm High Performance Back-illuminated Fiber Spectrometer

Our Product Introduction

for more products please visit us on spectralanalyser.com

Basic Information

- Place of Origin: CHINA
- Brand Name: JINSP
- Certification: ISO9001 CE
- Model Number: SR100B
- Minimum Order Quantity: 1
- Price: Negotiable
- Packaging Details: International Shipping Package
- Delivery Time: 30-400 working days
- Payment Terms: T/T, Western Union
- Supply Ability: 100PCS/30-40 days



Product Specification

- Spectral Range: 200nm - 1100nm
- Effective Pixels: 2048*64
- Pixel Size: 28.672*0.896mm
- Integration Time: 4ms~900s
- Dimensions: 180mm*120mm*50mm
- Weight: 1.2kg



Product Description

Backlit array CCD UV enhanced High Sensitivity Modular Spectrometer

The SR100B spectrometer stands out with its state-of-the-art back-illuminated CCD technology, providing unparalleled performance in the field. This versatile instrument is perfect for a variety of applications, from industrial processes to scientific research. The Hamamatsu S10420 CCD chip ensures a large sensing area and stable spectral readings, while the 2048x64 pixel format and 14x14μm pixel size guarantee remarkable quantum efficiency and sensitivity in the 200-1100nm range.

Customization options for gratings and slit widths are available to fine-tune spectral resolution and sensitivity. The SR100B's high-resolution optical path and sophisticated FPGA circuitry deliver noise-free, rapid signal processing, resulting in a highly reliable and stable spectrometer.

Product Parameters:

Detector	Chip Type	Back-illuminated cooling Hamamatsu S10420
	Effective Pixel	2048*64
	Pixel Size	14*14μm
	Sensing Area	28.672*0.896mm
Optical Parameters	Optical Design	F/4 cross type
	Numerical Aperture	0.13
	Focal Length	100mm
	Entrance Slit Width	10μm,25μm,50μm,100μm,200μm (customizable)
	Fiber Interface	SMA905,free space
Electrical Parameters	Integration Time	4ms~900s
	Data Output Interface	USB3.0,RS232,RS485,20pin connector
	ADC Bit Depth	16-bit
	Power Supply	5V
	Operating Current	<3.5A
Physical Parameters	Operating Temperature	10 ~40
	Storage Temperature	-20 ~60
	Operating Humidity	<90%RH(no condensation)
	Dimensions	180mm*120mm*50mm
	Weight	1.2kg

List of Product Models:

Model	Spectral Range(nm)	Resolution(nm)	Slit(μm)
SR100B-G21	200~1100	2.2	50
		1.5	25
		1.0	10
SR100B-G23	200~875	1.6	50
SR100B-G24	350~1025	1.0	25
		0.7	10
SR100B-G28	200~345	0.35	50
		0.2	25
		0.14	10
SR100B-G25	532~720(4900cm ⁻¹)*	13cm ⁻¹	50
SR100B-G26	638~830(3200cm ⁻¹)*	10cm ⁻¹	25
SR100B-G27	785~1080(3200cm ⁻¹)*	11cm ⁻¹	50

Note:The*are primarily designed for Raman applications, with the corresponding Raman.

Technical Characteristics:

High sensitivity-Fitted with area array back-illuminated detector with high quantum efficiency,optimizable ultraviolet band
 High resolution -Resolution<1.0nm@10μm (200~1100nm)
 High flexibility-180~1100nm,compatible with multiple interfaces including USB3.0,RS232 and RS485
 High reliability -Ultra-high SNR and excellent thermal

Typical Applications:


Detect absorption, transmittance and reflection Spectrum
 Light source and laser wavelength characterization
 OEM product module:Fluorescence spectrum,Raman spectrum,etc.

Technical Characteristics

- ✓ High sensitivity-Fitted with area array back-illuminated detector with high quantum efficiency, optimizable ultraviolet band
- ✓ High resolution -Resolution<1.0nm@10 μ m (200~1100nm)
- ✓ High flexibility-180~1100nm, compatible with multiple interfaces including USB3.0, RS232 and RS485
- ✓ High reliability -Ultra-high SNR and excellent thermal

Typical Applications

- ◆ Detect absorption, transmittance and reflection spectrum
- ◆ Light source and laser wavelength characterization
- ◆ OEM product module: Fluorescence spectrum, Raman spectrum, etc.

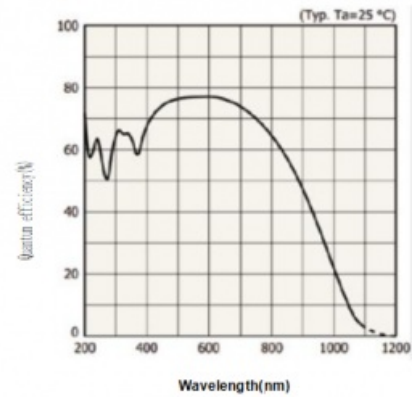

CCD Quantum Efficiency Curve

Company Introduction:

JINSP Company Limited, abbreviated as "JINSP", is a professional supplier with over 17 years of experience in spectral detection technology products, including Raman, FT-IR, LIBS technologies, etc. After 17 years of technology accumulation, the company's core key technologies have reached the international leading position at the level, and the cumulative number of patent applications exceeded 200.

In addition to its main headquarters located in the bustling city of Beijing, JINSP has established a fully owned subsidiary manufacturing facility situated in the province of Jiangsu, China.

JINSP Company received ISO9001:2015, ISO14001:2015, and ISO45001:2018 certifications. JINSP can provide required certifications, such as certification by the Ministry of Public Security or National Institute of Metrology, Environmental Level Certification, IP Level Certification, CE Certification, Transport Identification Report, EU ECAC certification, German ICT Security Testing, etc.

Company Profile



Exhibition



Certifications



FAQ

Q1: This is the first time I use, is it easy to operate?

A1:We will send you manual and guide video in English,it can teach you how to operate the spectrometer.Also our technicians will offer professional technical operation meetings.

Q2:Can you offer a operation training?

A2:Your technicians can come to our factory for a training. Jinsp engineers can go to your place for local support (installation , training, debugging, maintenance).

Q3:What's your website?

A3:You can visit:www.jinsptech.com

Q4:What about your quality assurance?

A4:We have a quality inspection team. All goods will go through quality inspection before shipment. We can send you pictures for inspection.



JINSP

JINSP Company Ltd.



8618620854039



phoebeyu@jinsptech.com



spectralanalyser.com

21st Floor, Building D, Tsinghua Tongfang Science and Technology Plaza, Haidian District, Beijing China